

A Novel Multiplex Tetra-Primer ARMS-PCR for the Simultaneous Genotyping of Six Single Nucleotide Polymorphisms Associated with Female Cancers

Chen Zhang^{1,9}, Ying Liu^{1,2,9}, Brian Z. Ring³, Kai Nie¹, Mengjie Yang¹, Miao Wang¹, Hongwei Shen¹, Xiyang Wu^{2*}, Xuejun Ma^{1*}

1 Key Laboratory for Medical Virology, Ministry of Health-National Institute for Viral Disease Control and Prevention, Chinese Center for Disease Control and Prevention, Beijing, China, 2 Department of Food Science and Engineering, College of Science and Engineering, Jinan University, Guangzhou, China, 3 Institute of Genomic and Personalized Medicine, School of Life Science and Technology, Huazhong University of Science and Technology, Wuhan, China

Abstract

Background: The tetra-primer amplification refractory mutation system PCR (T-ARMS-PCR) is a fast and economical means of assaying SNP's, requiring only PCR amplification and subsequent electrophoresis for the determination of genotypes. To improve the throughput and efficiency of T-ARMS-PCR, we combined T-ARMS-PCR with a chimeric primer-based temperature switch PCR (TSP) strategy, and used capillary electrophoresis (CE) for amplicon separation and identification. We assessed this process in the simultaneous genotyping of four breast cancer–and two cervical cancer risk–related SNPs.

Methods: A total of 24 T-ARMS-PCR primers, each 5'-tagged with a universal sequence and a pair of universal primers, were pooled together to amplify the 12 target alleles of 6 SNPs in 186 control female blood samples. Direct sequencing of all samples was also performed to assess the accuracy of this method.

Results: Of the 186 samples, as many as 11 amplicons can be produced in one single PCR and separated by CE. Genotyping results of the multiplex T-ARMS-PCR were in complete agreement with direct sequencing of all samples.

Conclusions: This novel multiplex T-ARMS-PCR method is the first reported method allowing one to genotype six SNPs in a single reaction with no post-PCR treatment other than electrophoresis. This method is reliable, fast, and easy to perform.

Citation: Zhang C, Liu Y, Ring BZ, Nie K, Yang M, et al. (2013) A Novel Multiplex Tetra-Primer ARMS-PCR for the Simultaneous Genotyping of Six Single Nucleotide Polymorphisms Associated with Female Cancers. PLoS ONE 8(4): e62126. doi:10.1371/journal.pone.0062126

Editor: Gayle E. Woloschak, Northwestern University Feinberg School of Medicine, United States of America

Received July 5, 2012; Accepted March 19, 2013; Published April 17, 2013

Copyright: © 2013 Zhang et al. This is an open-access article distributed under the terms of the Creative Commons Attribution License, which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited.

Funding: This work was supported by the China Mega-Project for Infectious Disease (2011ZX10004-001, 2012ZX10004-215 and 2013ZX10004-202). The funders had no role in study design, data collection and analysis, decision to publish, or preparation of the manuscript.

1

Competing Interests: The authors have declared that no competing interests exist.

- * E-mail: maxj2004@yahoo.com.cn (XM); tkentwu@jnu.edu.cn (XW)
- These authors contributed equally to this work.

Introduction

The role of single nucleotide polymorphisms (SNP's) in contributing to the variability between individuals in susceptibility to cancer[1], tumor growth and metastasis rate[2-4], as well as in treatment efficacy and adverse drug responses, has been well recognized[5,6]. Among the many methods that have been developed to genotype SNPs, the tetra-primer amplification refractory mutation system PCR (T-ARMS-PCR) has proved to be rapid, simple and economical [7–11]. Through combination of two outer primers and two allele-specific inner primers, genotyping requires only a single PCR followed by electrophoresis separation[8]. Multiplex PCR was incorporated in T-ARMS-PCR, using eight primers in one PCR, and is able to simultaneously detect two mutations [9]. Separately, chimericprimer-based multiplex PCR, which adds a universal 5' tag to the sequence specific primers for multiple targets, has been reported to improve the throughput and efficiency of the polymerase chain reaction[12]. With its high efficiency in detecting tens of different

PCR products in one reaction, the use of chimeric primer PCR has frequently been reported for use in mRNA quantification [13–15] and pathogen detection[16,17].

Breast cancer and cervical cancer have become the most frequently diagnosed cancers and the leading causes of cancer death among females[18]. Recent studies show that somatic variants in susceptibility regions are associated with the likelihood of occurrence of breast and gynecologic cancers [19-23]. SNPs were chosen for this study on the basis of reported associations with these cancers and having a reasonably high prevalence in Asian populations. Four low-penetrance variants for the prediction of risk of breast cancer were selected. SNPs rs4784227 [24] and rs3803662 [25] are located in the transcription factor TOX3; rs1219648 lies within FGFR2, which contributes to cell growth, invasiveness, motility, and angiogenesis [26]; rs889312 [27] is within MAP3K1, which is linked to cellular response to mitogens. Two variants associated with risk of cervical or ovarian cancer were selected. SNP rs750749 [21] is a polymorphisms in CD83, which is involved in immune recognition and antigen presentation; rs749292 in CYP19A1, which plays a key role in estrogen biosynthesis [22,23].

In this paper, we describe a novel multiplex T-ARMS-PCR allowing for the simultaneous genotyping of 6 SNPs (rs4784227, rs3803662, rs1219648, rs889312, rs750749 and rs749292) associated with breast and gynecologic cancers in a single tube using 24 chimeric primers and a pair of universal primers The use of chimeric primers and a temperature switch PCR (TSP) strategy were combined with T-ARMS-PCR to optimize the amplification parameters and improve the throughput of SNP genotyping. The combination of these different genotyping techniques demonstrates for the first time the ability of tetra-primer ARMS-PCR to reliably and efficiently detect six SNPs in a single reaction. Since more than 10 PCR products with different lengths need to be identified, capillary electrophoresis (CE) is used instead of agarose gel electrophoresis.

Materials and Methods

A total of 186 blood specimens from healthy Chinese female volunteers who were being monitored for potential hypertension were collected at community health centers in Wuhan, China during 2011 for this study. All aspects of the study were performed in accordance with the national ethics regulations and approved by the Institutional Review Boards of the Centre for Disease Control and Prevention 70 of China, as well as the Ethics Committee of Huazhong University of Science and Technology. Participants received "Written Informed Consent" of the study's purpose and of their right to keep information confidential. Written consent was obtained from all participants or their guardians.

Genomic DNA was extracted from 0.2 ml fresh peripheral blood samples by use of the Wizard Genomic DNA Purification Kit (Promega) according to manufacturer's instructions. Extracted DNA samples had a final concentration ranging from 55–365 ng/ μL .

To overcome the limitations of standard multiplex T ARMS-PCR methods, the proposed method was optimized in terms of primer design, PCR cycling conditions and in the utilization of chimeric primers and TSP strategy, as described in our previous reports in the detection of Influenza viruses and human hand foot and mouth associated pathogens[16,28]. A total of 24 chimeric primers each consisting of a gene-specific sequence with a universal tag sequence at the 5' end were used. The gene-specific portions of the primers were designed according to the requirements of T-ARMS-PCR. The specificity of allele-specific primers is conferred by the identity of the terminal 3' nucleotide with either the wild-type or the mutant allele, specificity is increased by the introduction of a deliberate mismatch at position -1 from the 3'-terminus. A pair of universal primers and six sets of T-ARMS-PCR chimeric primers were used for amplification. Detailed primer sequences and working concentrations for each SNP are listed in Table 1.

Genotyping of the assayed polymorphisms was performed by multiplex PCR amplification and fragment analysis. Six sets of T-ARMS-PCR primers for the amplification of twelve fragments of different sizes were pooled in a single 20 μl of reaction volume, which also contained 10 μL master mix of QIAgen Multiplex PCR kit, 50–100 ng of genomic DNA, and optimized concentrations of each primer (see Table 1). Multiplex PCR was performed using Bioer LifePro Thermal Cycler. An optimized temperature switch PCR (TSP) protocol, which uses four different annealing temperature was performed as follows: initial denaturation step of 95°C for 10 min, 3 cycles of 95°C for 30 s, 60°C for 30 s, and

 72° C for 45 s, 10 cycles of 95° C for 30 s, 58° C for 30 s, and 72° C for 45 s, 20 cycles of 95° C for 30 s, 68° C for 30 s, and 72° C for 45 s, 15 cycles of 95° C for 30 s, 55° C for 30 s, and 72° C for 45 s, followed by a final extension cycle at 72° C for 10 min, and then cooled to 4° C.

The multiplex PCR products were separated by QIAxcel DNA high-resolution gel cartridge (Qiagen) on QIAxcel system (Qiagen). DNA Size Marker of 25–450 bp (Qiagen) and Alignment Marker 15 bp/500 bp (Qiagen) were used in each QIAxcel runs and the size of the products was determined using the ScreenGel software (Qiagen). Because each of the amplicons was of a different length, the alleles were detected on the basis of the patterns of peak sizes.

A total of 186 samples were also sequenced in parallel with an ABI 3130 Genetic Analyzer (Applied Biosystems, USA) according to the BigDye Termination version 3.1 protocol in Invitrogen Corporation (Shanghai, China) to confirm the multiplex T-ARMS-PCR results using the outer primers listed in table 1 for each SNP.

Results

A total of 186 samples were typed with a multiplex T-ARMS-PCR assay and also typed in parallel with direct sequencing to assess the accuracy and efficiency of the assay. All the PCR products were well resolved and sized by CE and ScreenGel, allowing the easy identification of different genotypes. Two of the 186 samples had eleven unique amplicons, meaning those patients have only one homogenous SNP among the six tested loci. Electropherogram and gel image of these two samples (Fig. 1) show that CE on QIAXEL can clearly separate as many as 11 fragments in one sample. The accuracy of multiple PCR analysis of one sample was confirmed by direct sequencing (Fig. 2). Fragment sizes determined by QIAXEL based CE are listed in Table 2. Read lengths were 2 to 10 bp larger than the expected ones but this did not interfere with allele determination. Heterozygotes and homozygotes were unambiguously assigned from the CE profiles. No cross-reaction was observed. The genotypes scored from the multiplex assay were in 100% accordance with direct sequencing.

The genotype distribution and allele frequencies of each SNP are listed in table 3. The allele reported to be associated with risk of cancer occurrence is highlighted. The observed frequency of genotypes in this study was overall similar to that measured by HapMap for a Han Chinese population (HCB). If the HapMap frequencies were used to predict expected genotype counts in this study, then a comparison of this population to the HapMap HCB population showed significant divergence for rs4784227 in Tox3 and rs749292 in CYP19A1, in which the risk and non-risk alleles, respectively, are present at a significantly higher proportion than previously reported for a Chinese population. It is likely there is diversity within the Han population that is not yet captured by the HapMap studies. A supplementary table (Table S1) is provided to show the exact sets of alleles for all 6 SNPs found in each individual sample for further reference.

Discussion

Methods allowing low-cost, fast and reliable SNP determination are attracting increasing interest in the age of personalized medicine. It is widely acknowledged that using information from multiple SNP genotypes provides a more accurate risk assessment than that predicted by a single risk allele[29], therefore methods capable of identifying multiple genotypes, such as MALDI-TOF mass spectrometry [30,31] and hybridization-based [32–34] or

 Table 1. Primers for the multiplex tetra-primer ARMS-PCR for SNPs associated with breast and gynecologic cancers.

SNP	Primer	Primer Sequence	Working concentration(nM)	Length (bp
rs4784227	Forward inner primer(T):	AGGTGACACTATAGAATAAAAAGTC- CCAATTTGTAGTGTTTGaT ^{a,b,c}	40	322
	rs4784227-476T-2A			
	Reverse inner primer(C):	GTACGACTCACTATAGGGAAT-GGGAGTATTTACATCACAATAAT ${ m g}{ m G}^{a,b,c}$	40	297
	rs4784227-527G-2G			
	Forward outer primer:	AGGTGACACTATAGAATACTGA- CCCCTTTAGACACGG ^a	6.5	
	rs4784227-268			
	Reverse outer primer:	GTACGACTCACTATAGGGAGGGC- TTCAACACAGTCAGTTC ^a	6.5	
	rs4784227-760			
rs1219648	Forward inner primer(G):	AGGTGACACTATAGAATACACGC- CTATTTTACTTGACACAaG ^{a.b.c}	45	269
	rs1219648-278G-2A			
	Reverse inner primer(A):	GTACGACTCACTATAGGG- AAGCCATGGCCATCCTTGAAGAaT ^{a,b,c}	40	199
	rs1219648-323T-2A			
	Forward outer primer:	AGGTGACACTATAGAATACACAATGGCGCAGAATTA ^a	6.5	
	rs1219648-162			
	Reverse outer primer:	${\sf GTACGACTCACTATAGGGACTGGACAGGTCATTGTGGTG}^a$	7	
	rs1219648-509			
s3803662	Forward inner primer(C):	AGGTGACACTATAGAATACTCTCCTT- AATGCCTCTATAGCTGTaC ^{a.b.c}	45	259
	rs3803662-275C-2A			
	Reverse inner primer(T):	GTACGACTCACTATAGGGACCAC- AGTTTTATTCTTCGCTAACA ^{ab.c}	80	357
	rs3803662-324A-2C			
	Forward outer primer:	AGGTGACACTATAGAATAGTCCTTGGCTGTTCTGTGA ^a	11.5	
	rs3803662-5			
	Reverse outer primer:	GTACGACTCACTATAGGGATCCCCAAGGAGACAAAGGTA ^a	7	
	rs3803662-496			
rs889312	Forward inner primer(A):	AGGTGACACTATAGAATATGCCCC- TGCTGGAGAAAGtA ^{a,b,c}	60	213
	rs889312-382A-2T			
	Reverse inner primer(C):	GTACGACTCACTATAGGGATGATTTGT- AGTCTCTTAATTTGCACAaG ^{a.b.c}	75	340
	rs889312-428G-2A			
	Forward outer primer:	AGGTGACACTATAGAATAC- CTGGGTCCTTAGCATTCC ^a	11	
	rs889312-126			
	Reverse outer primer:	GTACGACTCACTATAGGGAATCT- GTGCCCTGAAGTGAGTAG ^a	8.5	
	rs889312-557			
s750749	Forward inner primer(C):	AGGTGACACTATAGAATAGAGGCTG- AGAACTAATAATATTTCC ^{a,b,c}	50	371
	rs750749-275C-2C			
	Reverse inner primer(T):	GTACGACTCACTATAGGGATCAG- ATGAAAACTGCTGTGTAAGaA ^{a,b,c}	75	224
	rs750749-325A-2A			
	Forward outer primer:	${\sf AGGTGACACTATAGAATACCCACAGAGTTGGAATGCTT}^a$	11	
	rs750749-139			
	Reverse outer primer:	GTACGACTCACTATAGGGATGCATAAAGTGGGTCCTGC ^a	7.5	

Table 1. Cont.

SNP	Primer	Primer Sequence	Working concentration(nM)	Length (bp)
rs749292	Forward inner primer(G):	${\sf AGGTGACACTATAGAATACCTTCTTCAAACCTCGGAGTaG}^{a,b,c}$	35	282
	rs749292-429G-2A			
	Reverse inner primer(A):	GTACGACTCACTATAGGGAGATAGA- AATTGTGCAGGAATCCaT ^{a.b.c}	25	330
	rs749292-472T-2A			
	Forward outer primer:	${\sf AGGTGACACTATAGAATATCCACGGACAGAGCAGG}^a$	4.5	
	rs749292-180			
	Reverse outer primer:	${\sf GTACGACTCACTATAGGGAGATTAGGGGACCTCCGTGA}^a$	6	
	rs749292-673			
Tag-F		AGGTGACACTATAGAATA ^d	500	
Tag-R		GTACGACTCACTATAGGGA ^d	500	

^aUnderlined is the universal tag sequence at the 5' end of the chimeric primers.

doi:10.1371/journal.pone.0062126.t001

enzyme-based[35] methods have been utilized. However these methods either require costly special equipment, such as mass spectrometer, or time consuming post-PCR operation.

Tetra-primer ARMS-PCR method has become one of the most commonly used methods for SNP genotyping. It requires only regular molecular biology equipment and eliminates the need of hybridization or additional enzymatic reactions. Although triplex and quadruplex PCR methods have been reported, there is limited use of multiplex T-ARMS-PCR in genotyping because of two key limitations. Firstly, the probability of finding primers with matched melting temperatures drastically drops when trying to combine the detection of several SNPs in a single reaction. Secondly, the

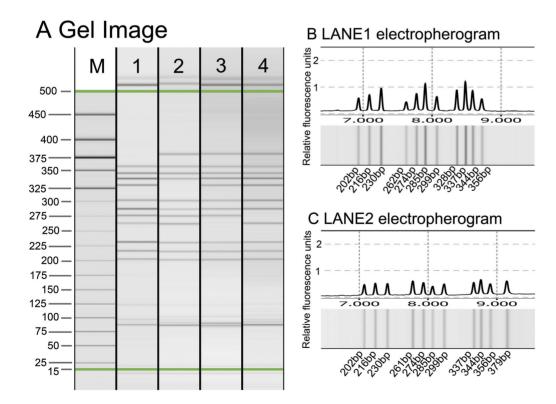


Figure 1. Electrophoresis result of 4 samples with most number of amplicons. A: gel image of the 4 samples. B. Gel image and electropherogram of sample in lane1. 11 bands ranging from 202 to 356 bp were observed, indicating 5 heterogenous and 1 homogenous SNPs were identified C: Gel image and electropherogram of sample in lane2. Another combination of 5 heterogenous and 1 homogenous SNPs were observed.

doi:10.1371/journal.pone.0062126.g001

^bAllele-specific nucleotides at the 3' end are indicated in bold letters.

^cSpecificity is increased by the introduction of a deliberate mismatch at position -1 of the polymorphism site, indicated by bold lower case letters.

 $^{^{\}mathrm{d}}\mathrm{A}$ pair of universal primers annealing to the 5' portion of each chimeric primer.

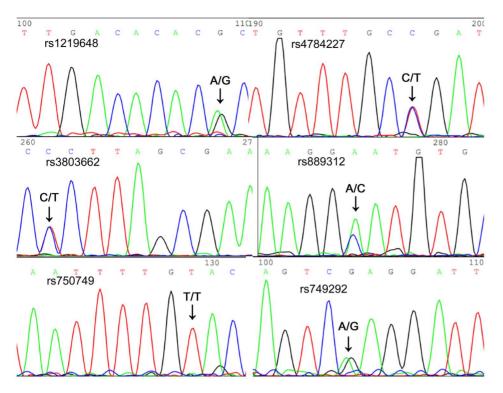


Figure 2. Direct sequencing result of sample #1101-291(corresponding to the sample in lane 1, Fig.1B), which showed the combination of 5 heterogenous and 1 homogenous SNPs. doi:10.1371/journal.pone.0062126.q002

resulting pool of amplicons requires sufficient length intervals between neighboring bands in electrophoresis to facilitate separation. By combining T-ARMS-PCR with a chimeric primer-based temperature switch PCR strategy we largely sidestep these limitations. Our method demonstrates for the first time the ability of tetra-primer ARMS-PCR to easily detect six SNPs in a single reaction.

The reliability of the method was illustrated by typing 186 clinical blood samples in parallel with direct sequencing, and a 100% consistency between the two methods was obtained. This

proof-of-concept study thus establishes a rapid, reproducible, and cost effective method for the detection of multiplex SNPs, since CE by QIAXCEL is capable of resolving amplicons with as little as 5 bp size difference, the smallest size difference in this test was 10 bp. As the read lengths in this assay has a standard deviation of 0.8 to 1.3 (Table 2), determination of allele is thus not interfered.

The use of chimeric primers and biphasic temperature switch in the annealing process decreases the difference in amplification efficiency among amplicons. During the first few PCR cycles, amplification is carried out by allele-specific chimeric primers. In

Table 2. Analyses of read lengths in electrophoresis.

		Min. read	Max. read	Average	Sd.of read	
	Theoretical length	length	length	read length	length	99% CI of mean
rs1219648-A(199)	199	201	206	203.09	1.18	(200–206)
rs889312-A(213)	213	215	220	216.48	1.17	(213–219)
rs750749-T(224)	224	229	234	231.21	1.13	(228–234)
rs3803662-C(259)	259	259	265	261.98	1.09	(259–265)
rs1219648-G(269)	269	273	277	274.74	1	(272–277)
rs749292-G(282)	282	282	289	285.91	1.12	(283–289)
rs4784227-C(297)	297	298	305	300.16	1.02	(298–303)
rs4784227-T(322)	322	327	330	328.39	0.8	(326–330)
rs749292-A(330)	330	336	339	337.5	0.63	(336–339)
rs889312-C(340)	340	342	347	344.73	0.82	(343–347)
rs3803662-T(357)	357	352	360	356.1308	1.34	(353–360)
rs750749-C(371)	371	377	382	379.3372	0.86	(377-382)

doi:10.1371/journal.pone.0062126.t002

Table 3. Genotypes and allele frequencies of tested samples in this study.

	genotype	N (n=186)	Frequency (%)	Expected genotype frequency under HWE(%)	HapMap frequency HCB (%)	Allele	Allele frequency (%) HW ($ ho$ -value) ^{b}	НW (<i>p</i> -value) ^b
rs4784227	CC	93	50.00	51.52	66.7	O	77.17	0.31
	CJ	81	43.55	40.52	24.4	-	28.23	
	F	12	6.45	7.97	8.9			
rs1219648	AA	99	35.48	34.66	40	A	58.87	0.64
	AG	87	46.77	48.43	46.7	ט	41.13	
	99	33	17.74	16.92	13.3			
rs3803662	SS	13	6.99	8.59	6.7	O	29.30	0.29
	CJ	83	44.62	41.43	44.4	-	70.70	
	F	06	48.39	49.98	48.9			
rs889312	CC	41	22.04	22.64	27.3	U	47.58	0.74
	AC	95	51.08	49.88	45.5	۷	52.42	
	AA	50	26.88	27.48	27.3			
rs750749	F	143	76.88	77.74	80	-	88.17	0.26
	C	42	22.58	20.86	17.8	O	11.83	
	SS	-	0.54	1.40	2.2			
rs749292	99	62	33.33	30.07	8.9	9	54.84	0.07
	GA	80	43.01	49.53	55.6	4	45.16	
	AA	44	23.66	20.40	35.6			

^aincreased risk allele in bold ^bp value was calculated using chi test doi:10.1371/journal.pone.0062126.t003 later stages of PCR, amplification is predominantly carried out by universal primers, so that all targets in this multiplex PCR system are amplified in an unbiased manner by a single pair of universal primers. This reduces the occurrence of biased and partial amplification, minimizes non-specific reactions, and reduces the need for optimization of each individual PCR assay.

To assess the errors resulting from the utilization of multi-band electrophoregrams to distinguish amplicons, the read length of each band was compared with the theoretical lengths calculated from primer-alignment (Table 2). No overlapping in read length range between any two of the amplicons was observed, additionally, there is no overlap of the 99% confidence intervals of the observed average read length. It is notable that band intensities can be subject to several factors including genomic DNA quality and PCR reagents quality; we have found that 50-100 ng/µL DNA is optimal to provide sufficiently clear and bright bands with minimal background (data not shown). In addition, unlike most other reported T-ARMS-PCR methods, a deliberated mismatch at position -1 from the 3' terminus was incorporated into both inner primers and was sufficiently specific for the differential detection of two alleles for each SNP. Due to the limitation of size discrimination between PCR products, PCR primer design may be restricted to some extent, making multiple T-ARMS-PCR difficult to type those SNPs that are located closer than 20 bp to each other.

Two distinct advantages of the multiple ARMS-PCR method are the short assay time and the low costs, even for assaying large

References

- Erichsen HC, Chanock SJ (2004) SNPs in cancer research and treatment. Br J Cancer 90: 747–751.
- Cully M, You H, Levine AJ, Mak TW (2006) Beyond PTEN mutations: the PI3K pathway as an integrator of multiple inputs during tumorigenesis. Nat Rev Cancer 6: 184–192.
- Hollstein M, Sidransky D, Vogelstein B, Harris CC (1991) p53 mutations in human cancers. Science 253: 49–53.
- 4. Hanahan D, Weinberg RA (2000) The hallmarks of cancer. Cell 100: 57–70.
- Sjoqvist F (1999) The past, present and future of clinical pharmacology. Eur J Clin Pharmacol 55: 553–557.
- Wiechec E, Hansen LL (2009) The effect of genetic variability on drug response in conventional breast cancer treatment. Eur J Pharmacol 625: 122–130.
- Newton CR, Graham A, Heptinstall LE, Powell SJ, Summers C, et al. (1989) Analysis of any point mutation in DNA. The amplification refractory mutation system (ARMS). Nucleic Acids Res 17: 2503–2516.
- Ye S, Humphries S, Green F (1992) Allele specific amplification by tetra-primer PCR. Nucleic Acids Res 20: 1152.
- Piccioli P, Serra M, Gismondi V, Pedemonte S, Loiacono F, et al. (2006) Multiplex tetra-primer amplification refractory mutation system PCR to detect 6 common germline mutations of the MUTYH gene associated with polyposis and colorectal cancer. Clin Chem 52: 739–743.
- Balbi G, Ferrera F, Rizzi M, Piccioli P, Morabito A, et al. (2007) Association of 318 C/T and +49 A/G cytotoxic T lymphocyte antigen-4 (CTLA-4) gene polymorphisms with a clinical subset of Italian patients with systemic sclerosis. Clin Exp Immunol 149: 40–47.
- Peruzzi B, Serra M, Pescucci C, Sica M, Lastraioli S, et al. (2010) Easy genotyping of complement C3 'slow' and 'fast' allotypes by tetra-primer amplification refractory mutation system PCR. Mol Cell Probes 24: 401–402.
- Shuber AP, Grondin VJ, Klinger KW (1995) A simplified procedure for developing multiplex PCRs. Genome Res 5: 488–493.
- Wang QX, Li K, Zhou YX, Xiao JH (2009) [Development, optimization and application of the expression analysis platform based on multiplex quantitative RT-PCR using fluorescent universal primers]. Yi Chuan 31: 552–561.
- Beck O, Seidl C, Lehrnbecher T, Kreyenberg H, Schwabe D, et al. (2006) Quantification of chimerism within peripheral blood, bone marrow and purified leukocyte subsets: comparison of singleplex and multiplex PCR amplification of short tandem repeat (STR) loci. Eur J Haematol 76: 237–244.
- Dupont M, Goldsborough A, Levayer T, Savare J, Rey JM, et al. (2002) Multiplex fluorescent RT-PCR to quantify leukemic fusion transcripts. Biotechniques. 33: 158–160, 162, 164
- Hu X, Zhang Y, Zhou X, Xu B, Yang M, et al. (2012) Simultaneously typing nine serotypes of enteroviruses associated with hand, foot, and mouth disease by a GeXP analyzer-based multiplex reverse transcription-PCR assay. J Clin Microbiol 50: 288–293.

numbers of specimens. The proposed method, only involving conventional PCR with a CE, can be performed within 3.5 hours with minimal hands-on effort. After extraction of genomic DNA, the subsequent steps can be completed in a single reaction tube, allowing for the ready analysis of multiple samples in a single run for high-throughput screening. This assay consumes only standard PCR reagents and electrophoresis cartridges; the costs in this study were only 2 US\$ for the simultaneous detection of six SNPs per sample.

To our knowledge, the proposed method is the first to detect six SNPs in a single reaction using tetra-primer ARMS-PCR. The novel multiplex tetra-primer ARMS-PCR method developed in this study has significant potential to be widely applicable in both commercial and clinical settings for the screening of multiple SNPs.

Supporting Information

Table S1 Alleles of 6 SNPs of each individual sample in this study.

(DOCX)

Author Contributions

Sample collection and DNA extraction: MW HS. Conceived and designed the experiments: CZ XM XW. Performed the experiments: YL CZ. Analyzed the data: YL CZ KN. Contributed reagents/materials/analysis tools: YL CZ BR MY. Wrote the paper: YL CZ BR XM.

- Hofmann MA (2003) Construction of an infectious chimeric classical swine fever virus containing the 5'UTR of bovine viral diarrhea virus, and its application as a universal internal positive control in real-time RT-PCR. J Virol Methods 114: 77–90.
- Jemal A, Bray F, Center MM, Ferlay J, Ward E, et al. (2011) Global cancer statistics. CA Cancer J Clin 61: 69–90.
- Easton DF, Pooley KA, Dunning AM, Pharoah PD, Thompson D, et al. (2007) Genome-wide association study identifies novel breast cancer susceptibility loci. Nature 447: 1087–1093.
- Garcia-Closas M, Chanock S (2008) Genetic susceptibility loci for breast cancer by estrogen receptor status. Clin Cancer Res 14: 8000–8009.
- Yu KJ, Rader JS, Borecki I, Zhang Z, Hildesheim A (2009) CD83 polymorphisms and cervical cancer risk. Gynecol Oncol 114: 319–322.
- Setiawan VW, Doherty JA, Shu XO, Akbari MR, Chen C, et al. (2009) Two
 estrogen-related variants in CYP19A1 and endometrial cancer risk: a pooled
 analysis in the Epidemiology of Endometrial Cancer Consortium. Cancer
 Epidemiol Biomarkers Prev 18: 242–247.
- Goodman MT, Lurie G, Thompson PJ, McDuffie KE, Carney ME (2008)
 Association of two common single-nucleotide polymorphisms in the CYP19A1 locus and ovarian cancer risk. Endocr Relat Cancer 15: 1055–1060.
- Long J, Cai Q, Shu XO, Qu S, Li C, et al. (2010) Identification of a functional genetic variant at 16q12.1 for breast cancer risk: results from the Asia Breast Cancer Consortium. PLoS Genet 6: e1001002.
- Reeves GK, Travis RC, Green J, Bull D, Tipper S, et al. (2010) Incidence of breast cancer and its subtypes in relation to individual and multiple lowpenetrance genetic susceptibility loci. JAMA 304: 426–434.
- Jia C, Cai Y, Ma Y, Fu D (2010) Quantitative assessment of the effect of FGFR2 gene polymorphism on the risk of breast cancer. Breast Cancer Res Treat 124: 521–528.
- Lu PH, Yang J, Li C, Wei MX, Shen W, et al. (2011) Association between mitogen-activated protein kinase kinase kinase 1 rs889312 polymorphism and breast cancer risk: evidence from 59,977 subjects. Breast Cancer Res Treat 126: 663–670.
- Qin M, Wang DY, Huang F, Nie K, Qu M, et al. (2010) Detection of pandemic influenza A H1N1 virus by multiplex reverse transcription-PCR with a GeXP analyzer. J Virol Methods 168: 255–258.
- Harlid S, Ivarsson MI, Butt S, Grzybowska E, Eyfjord JE, et al. (2012) Combined effect of low-penetrant SNPs on breast cancer risk. Br J Cancer 106: 389–396.
- Millis MP (2011) Medium-throughput SNP genotyping using mass spectrometry: multiplex SNP genotyping using the iPLEX(R) Gold assay. Methods Mol Biol 700: 61–76
- 31. Xiu-Cheng Fan A, Garritsen HS, Tarhouny SE, Morris M, Hahn S, et al. (2008) A rapid and accurate approach to identify single nucleotide polymorphisms of

- mitochondrial DNA using MALDI-TOF mass spectrometry. Clin Chem Lab Med 46: 299-305.
- Shi X, Tang C, Zhou D, Sun J, Lu Z (2009) PCR-product microarray based on polyacrylic acid-modified surface for SNP genotyping. Electrophoresis 30: 1286– 1296
- Yanagawa T, Koga H (2009) Modified multiple primer extension method. Methods Mol Biol 578: 425–435.
- 34. Tebbutt SJ, Ruan J (2008) Combining multiple PCR primer pairs for each amplicon can improve SNP genotyping accuracy by reducing allelic drop-out. Biotechniques. 45: 637–638, 640, 642 passim
- Toubanaki DK, Christopoulos TK, Ioannou PC, Flordellis CS (2009) Identification of single-nucleotide polymorphisms by the oligonucleotide ligation reaction: a DNA biosensor for simultaneous visual detection of both alleles. Anal Chem 81: 218–224.